





# bio needle<sup>®</sup> H24

Adjustable needle size

## Needle Depth 0-1.5mm

The H24 allows users to choose appropriate needle depth depending on their skin type and treating areas to address different skincare concerns.



Forehead



Cheek



Eye



Body



## Hygiene and dual protection


The plastic shell can assure the needles and the glass ampoule are safe and clean.

The glass ampoule is detachable which could be replaced with standard glass ampoule on the market.



## Evenly & Efficiently

While the H24 is working on the skin surface, the skin serum will discharge evenly from 24 tiny holes and central axis, follow the needles and reach different skin layers. The H24 plays well with different skin serums to improve the skin problem

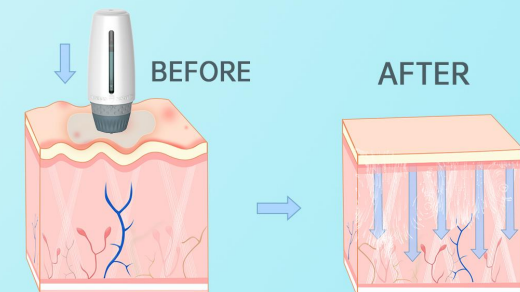

## Hygiene and dual protection

The plastic shell can assure the needles and the glass ampoule are safe and clean.

The glass ampoule is detachable which could be replaced with standard glass ampoule on the market.



1




2



3



4

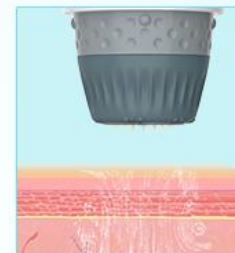




## Needle Depth 0-1.5mm

The H24 allows users to choose appropriate needle depth depending on their skin type and treating areas to address different skincare concerns.



## Evenly & Efficiently

While the H24 is working on the skin surface, the skin serum will discharge evenly from 24 tiny holes and central axis, follow the needles and reach different skin layers. The H24 plays well with different skin serums to improve the skin problem

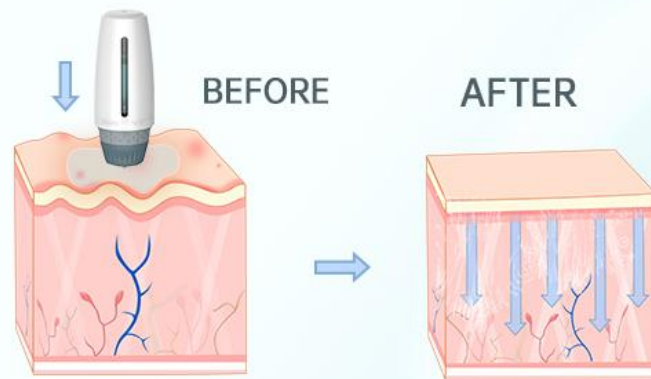

**Instructions:**

1. Remove the head from the case by turning the head counterclockwise.
2. Take off the head from the vial.
3. Put the essence of the liquid into the vial and then assemble the head back to the vial.
4. Lock the head with vial back to the case.
5. Take off the cap from the head.
6. Adjust the needle length by turning the scale ring.  
(The needle length will be the number pointing by triangle.)



OPEN





Product Parameter			
Name	H24	Capacity of Vial	10ML
Handle Material	PE+Grass	Instructions	Single use
Available Needle Length	0-1.5mm	Package Weight	50g
Package Size	101*67*33mm	Packing List	①Carton ②Tray ③H24