



# iBEAUTY Mate

Intelligent, multifunctional device  
for cosmetic procedures




## Welcome to the world of **iBeauty Mate**

the ultimate sidekick on your journey to beauty and self-confidence.  
Say goodbye to costly cosmetic procedures and say hello to a revolutionary  
all-in-one device that brings professional salon treatments right to your  
fingertips.



**Effective and safe  
treatments with IBEAUTY MATE**

IBEAUTY MATE allows both the novice and a professional cosmetologist to carry out the following salon-quality, safe and effective treatments in the comfort of their own home.

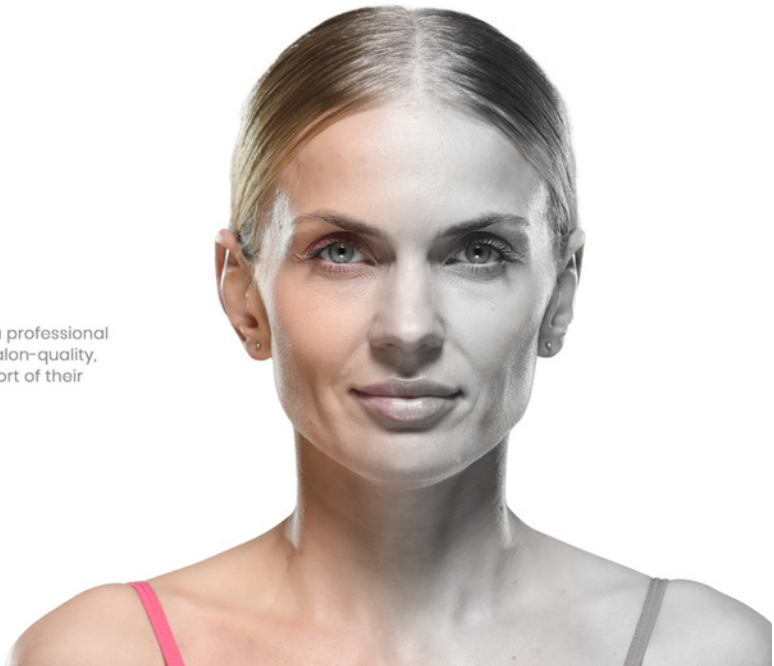

# The IBEAUTY MATE

device consists of the following components



- **A spray gun unit designed for skin treatments such as oxygen infusion and oxyabrasion** – an excellent way to cleanse the skin, remove dead skin cells and prepare the tissue for active ingredients intake.

#### Technical data

Pump capacity	25-28 l/min
Pump pressure	0-12 bar
Output maximum pressure	2 bar
Head dimensions	15 x 7 x 3 cm
Head weight	110 g
Tank capacity	7 cm <sup>3</sup>
Head material	aluminium
Hose length	1,8 m



- **A nano-needle head unit intended for needle free mesotherapy treatments** – an effective way of slowing down the progressive ageing process and eliminating wrinkles, discolouration, acne scars, under-eye dark circles, skin laxity, dryness, and lack of skin elasticity, lipodystrophy, stretch marks or excessive hair loss.

#### Technical data

Cartridge type	with 0,1 mm nano-needles level 1: 25-30 insertions/sec level 2: 35-40 insertions/sec
Frequency	level 3: 45-50 insertions/sec level 4: 55-60 insertions/sec level 5: 65-70 insertions/sec
Line feed depth	0 - 3 mm
Head dimensions	15 x 2 cm
Head weight	77 g
Head material	aluminium
Cable length	1,8 m



- **A head unit for diamond microdermabrasion treatment** – a treatment to instantly smooth, brighten and cleanse the skin, improve its tone and provide anti-ageing benefits.
- **A set of abrasive head units designed for microdermabrasion treatments** – treatments to effectively improve skin condition and eliminate imperfections such as wrinkles, discolouration, scars or acne, blackheads and unaesthetic enlarged pores.

Technical data	
Pump capacity	12-15 l/min
Pump pressure	0-2,2 bar
Types of abrasive tips	7 mm (D-240)
	8 mm (D-200)
	14 mm (D-180)
	14 mm (D-180)
	14 mm (D-120)
	17 mm (D-240)
	17 mm (D-180)
	20 mm (D-180)
	20 mm (D-75) FOR FEET ONLY
	22 mm (D-120)
Head dimensions	1,5 x 11 cm
Filtration type	pre-filter, liquid filter on the hose
Hose length	1,8 m



- **Two various-sized heads for endermic massage** – ideal for smaller and larger body areas.

#### Technical data

Pump capacity	12-15 l/min
Pump pressure	0-2,2 bar
Under-pressure	0- 53 cmHg
Large head diameter	12 cm
Large head weight	700 g
Small head diameter	7,5 cm
Small head weight	280 g
Head material	Poliacetal
Roller material	Stainless steel
Filtration type	pre-filter + liquid filter on the hose
Hose length	1,8 m





- **A head unit designed for superficial kinetic massage treatments** – treatments with analgesic and anti-inflammatory effects that stimulates blood circulation, eliminates cellulite and relaxes.

#### Technical data

Pump capacity	25-28 l/min
Output maximum pressure	2 bar
Pump pressure	0-12 bar
Head dimensions	3,7 x 19 cm
Head weight	650 g
Head material	Poliacetal, stainless steel, aluminium
Filtration type	condenser on the hose
Hose length	1,8 m



- **A head unit designed for LED light therapy treatments** – treatments to accelerate skin regeneration processes and alleviate erythema, scar irritation and hyperpigmentation.

#### Technical data

Number of diodes	7
Current for a single diode	18 mA
Colours	white, blue, green, red, orange, yellow, purple
Dimensions	3,7 x 19 cm
Head weight	130 g
Head material	Poliaceta
Cable length	1,8 m



## How does **iBeauty Mate WORK**

The innovative features of the IBEAUTY MATE multifunctional cosmetic device are encapsulated in the following functions.



## Personalised user profiles

IBEAUTY MATE allows you to create several user profiles – you can use it alone, or with all your family or friends, and get the results you desire according to your requirements.

What's more, thanks to an intelligent matrix, once you have indicated your contraindications, the device independently selects safe and effective treatments for you in relation to them and prevents you from performing inadvisable procedures.

You can also save your favourite treatments in your user profiles for convenience and get information on when and how you should perform the subsequent procedure.



## Professional step-by-step treatments for various body parts

The IBEAUTY MATE software suggests treatments for specific body zones – from the scalp and face, through the neck, décolleté and body, to the buttocks, legs, arms, hands, and feet. Therefore, you can focus your treatment on specific body areas and maximise your desired results.




Once you have selected a specific treatment for a particular body part, you see detailed information about the treatments, technology, and any accessories you may require to perform the treatment on the IBEAUTY MATE screen.

The IBEAUTY MATE device, using easy-to-understand graphics, demonstrates how to connect a particular head and use it. Then all that is left to do is to familiarise yourself with the childlike step-by-step instructions in text and animated form.

The IBEAUTY MATE device also proposes average parameters, doses, and timing to ensure even more safety during use. Finally, IBEAUTY MATE recommends actions for you to maximise the effects of the treatment. And that's it! You can enjoy great results – just as you would after treatment with a professional beautician.



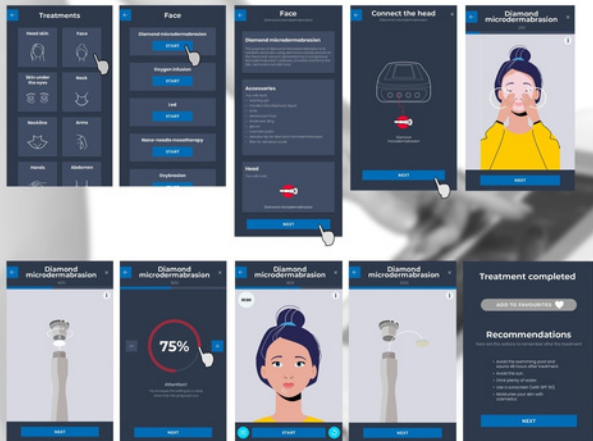

## Personalised user profiles

Besides treatments for specific parts of the body, IBEAUTY MATE also offers you perfectly composed, in terms of parameters, timing and frequency, dedicated treatment programmes designed to combat specific aesthetic problems, such as complete therapies aimed at eliminating stretch marks, fighting cellulite, treatments to accelerate hair growth or wrinkle therapy.

## BEAUTY MATE PRO –dedicated to professionals

The parameters of IBEAUTY MATE are the equivalent of professional equipment, so it can successfully substitute similar machinery found in beauty salons. With this application in mind, we have just created the IBEAUTY MATE PRO programme, geared towards time-saving and dedicated exclusively to professionals.

IBEAUTY MATE PRO is the ideal software for professional beauticians working in salons or travelling to clients' homes. This part of the IBEAUTY MATE device is devoid of instructions and any unnecessary information for the experienced professional and allows you to go straight to business. It enables professional cosmetologists, aestheticians, and cosmetic practitioners to save time and increase revenue without compromising quality and safety.





# Rejuvenate

both your skin and your  
life today!

IBEAUTY MATE gives them the ultimate go-to for obtaining professional results without paying for expert treatments.

It's the ideal accompaniment to your daily skincare routine, so give this solution a tryout and give yourself the time, money and confidence!



## Key benefits

IBEAUTY MATE is a safe and practical all-in-one device to bring you closer to getting the top beauty salon-like results without the hassle of finding the time and extra resources for regular treatments.

It's the perfect consensus between professional outcomes and savings in time and money!

So let us tell you about it in a little more detail.

What are the main advantages of the IBEAUTY MATE device?





## Convenience

IBEAUTY MATE is a compact and handy sidekick on the road to beauty, relaxation and confidence. We have ensured that the design and structure of our device allow it to travel with you and be easy to use.



## Saving time and money

Anyone who lives and works in a big city knows how difficult it is to juggle all your responsibilities while still remembering to regularly pay visits to the cosmetologist or beautician. As a long-standing supplier to such venues, we also recognise this issue. That is why, we have developed IBEAUTY MATE, a professional device that will save you time!



## Safety

We are a European manufacturer of professional cosmetic equipment with many years of experience. Therefore, we are confident in the safety and effectiveness of our revolutionary device.



## Effectiveness

As a manufacturer of professional devices for use in beauty salons and aesthetic practices, we boast a vast experience in developing professional devices that deliver actual results.

The IBEAUTY MATE is not another fleeting beauty world trend but a genuine, professional cosmetological device condensed to a compact size and optimally error-proofed to deliver actual results.

← Treatment

Head skin



Skin under the eyes



Neckline



IBEAUTY  
MATE

# Technical data

General	
Power consumption	150 W
Voltage	220-230 V
Frequency	50-60 Hz
Display	7" LCD, touch-sensitive IPS 1024 x 600 px
Ventilation method	Air fan
Base dimensions	38,5 x 21 x 21,5 (cm)
Weight	7,4 kg
Standards / Certification	CE

Nano-needle mesotherapy	
Cartridge type	with 0,1 mm nano-needles level 1: 25-30 insertions/sec level 2: 35-40 insertions/sec level 3: 45-50 insertions/sec level 4: 55-60 insertions/sec level 5: 65-70 insertions/sec
Frequency	
Line feed depth	0 - 3 mm
Head dimensions	15 x 2 cm
Head weight	77 g
Head material	aluminium
Cable length	1,8 m

Diamond microdermabrasion	
Pump capacity	12-15 l/min
Pump pressure	0-2,2 bar
	7 mm (D-240)
	8 mm (D-200)
	14 mm (D-180)
	14 mm (D-180)
	14 mm (D-120)
Types of abrasive tips	17 mm (D-240)
	17 mm (D-180)
	20 mm (D-180)
	20 mm (D-75) FOR FEET ONLY
	22 mm (D-120)
Head dimensions	1,5 x 11 cm
Filtration type	pre-filter, liquid filter on the hose
Hose length	1,8 m

Oxygen infusion	
Pump capacity	25-28 l/min
Pump pressure	0-12 bar
Output maximum pressure	2 bar
Head dimensions	15 x 7 x 3 cm
Head weight	110 g
Tank capacity	7 cm <sup>3</sup>
Head material	aluminium
Hose length	1,8 m

Superficial kinetic massage	
Pump capacity	25-28 l/min
Output maximum pressure	2 bar
Pump pressure	0-12 bar
Head dimensions	3,7 x 19 cm
Head weight	650 g
Head material	Polycetal, stainless steel, aluminium
Filtration type	condenser on the hose
Hose length	1,8 m

Endermossage	
Pump capacity	12-15 l/min
Pump pressure	0-2,2 bar
Under-pressure	0- 53 cmHg
Large head diameter	12 cm
Large head weight	700 g
Small head diameter	7,5 cm
Small head weight	280 g
Head material	Polycetal
Roller material	Stainless steel
Filtration type	pre-filter + liquid filter on the hose
Hose length	1,8 m

LED	
Number of diodes	7
Current for a single diode	18 mA
Colours	white, blue, green, red, orange, yellow, purple
Dimensions	3,7 x 19 cm
Head weight	130 g
Head material	Polycetal
Cable length	1,8 m

Oxybrasion	
Pump capacity	25-28 l/min
Pump pressure	0-12 bar
Output maximum pressure	2 bar
Head dimensions	15 x 7 x 3 cm
Head weight	110 g
Tank capacity	7 cm <sup>3</sup>
Head material	aluminium
Hose length	1,8 m



# Distributor in Greece

## Contact:

Sigma Medical Co - Medical.gr  
Ippokratous 142 - 19004 Spata  
Ag. Thoma 22 - 11527 Athens  
K. Melenikou 7 - 54635 Thessaloniki

Web: [medical.gr](http://medical.gr)  
email: [sales@medical.gr](mailto:sales@medical.gr)  
Tel. +30 210 6630520

iBEAUTY  
Mate