

### HEINE Large Battery Handle



Dimensions: 165 mm x 37 mm dia. Weight: 212g without batteries, 505g with batteries.

A large handle for 2 alkali-manganese batteries IEC LR20/D cells. Designed for countries where larger batteries are still preferred. (Only suitable for certain instrument cases.) Batteries not included.

<b>Large battery handle</b>	<b>2.5V</b>
for 2 alkali-manganese batteries, size D	<b>X-001.99.120</b>

### HEINE mini 3000® Battery Handle



Innovative, compact battery handle in a high-quality, modern design. Suitable for all mini3000 instruments.

- ∴ **Oval, compact, stylish shape.** Lies ergonomically in the hand.
- ∴ **High-quality handle: Chrome-finish upper section/refined plastic.** Shockproof, sturdy, non-slip.
- ∴ **Attachment clip with integrated on/off switch.** Secure. Switches off automatically when returned to the pocket. 20.000 switch cycles guaranteed.
- ∴ **Screw-on end cap.** For easy battery change.
- ∴ **Threaded connector made of metal.** Prevents wear.
- ∴ **Replaceable batteries.** Size AA, or optional rechargeable battery with mini NT table charger.
- ∴ **Complete range of mini 3000 instruments.** Versatile range.
- ∴ **Registered Design.**
- ∴ **Rechargeable mini 3000 handle.** The mini 3000 battery handle can be retro-fitted with the mini rechargeable battery 2.5V NiMH and the bottom insert to update to the mini 3000 rechargeable handle. Fully-automatic charging with the mini NT table charger.

<b>mini 3000 Battery Handle</b>	<b>2.5V</b>
with two dry batteries IEC LR6 (AA size)	<b>D-001.79.021</b>



All non-rechargeable mini handles are fitted with batteries (IEC LR 6 AA size). mini3000 handles can be supplied in black or blue. Please specify the colour you prefer when ordering. If not otherwise specified, we will supply black instruments.



## HEINE® mini NT Charger Set for mini3000®

### Compact charging technology



∴ Fast charge capability

∴ Automatic charging

∴ Charge indicator


**mini NT Charger set 2.5V.** Consists of the mini NT table charger, 2 mini NT rechargeable batteries and 2 bottom inserts for the mini 3000 rechargeable handles. With this set, existing mini 3000 handles can be updated to the 2.5V rechargeable system. Optional adaptor set for mini 2000.  
mini NT table charger: Width 112mm, Height 50mm, Depth 52mm, Weight 270g.

- ∴ **Innovative, compact design.** Practical and space-saving.
- ∴ **Two independent charging wells.** Complete functionality for one or two handles at a time.
- ∴ **Fast charge feature.** Typ. 4 hours to recharge a mini3000 NiMH 2.5V rechargeable battery.
- ∴ **Automatic charging.** Independent of the state of charge of the rechargeable battery. Guarantees extended battery life.
- ∴ **Cannot be overcharged.** Charger switches off automatically when the battery is fully charged.
- ∴ **Charge status indicator.** Light pulses from the charging well to confirm the charging mode. Continuous light indicates that the charging process has been completed.
- ∴ **No memory effect.** Full power at any time without loss of capacity.


**Please note:** Suitable only for mini3000 2.5V NiMH rechargeable battery.



mini NT table charger and accessories	
mini NT Charger, only [01]	X-001.99.484
mini3000 Rechargeable handle complete* [02]	D-001.79.022
Rechargeable battery NiMH 2Z, only [03]	X-001.99.487
Bottom insert for mini 3000 Rechargeable Handle, only [04]	X-001.99.488

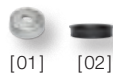

mini NT Adaptor Set for 1 mini 2000 Handle	
This set contains:	
1 bottom insert [01] and 1 adaptor for mini 2000 [02]. Together with a rechargeable battery it will allow the upgrade of one existing mini 2000 handle.	X-001.99.493

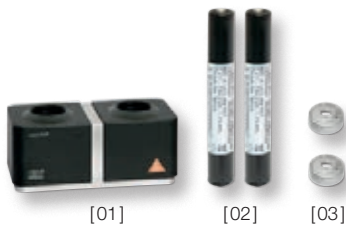

mini NT Charger Set for 2 mini 3000 Handles	
With this set, existing mini 3000 Handles can be updated to the 2.5V rechargeable system. It contains:	
mini NT table charger [01], 2 rechargeable batteries NiMH 2Z [02], 2 bottom inserts for mini 3000 rechargeable handle [03]	X-001.99.485



\*HEINE mini 3000 pocket instruments – a complete range of quality diagnostic instruments for your pocket. Available in black or blue. Please specify the colour you prefer when ordering. If not otherwise specified, we will supply black instruments.