

Digital Blood Pressure Monitor --sphygmomanometer Model : CONTEC 08A



Features:

1.Easy Use

With just the touch of one button, the monitor auto inflates for quick, easy measurements.

2.Colorful, high resolution , large LCD display panel shows readings of blood pressure, pulse rate and optional SPO2 displayed simultaneously.

Powerful functions

- 1.Three-user data storage,maximum 100 records for each user, covering the need of all your family***
- 2.Three mode for adult, pediatric and neonatal selectable***
- 3.You can set alarms so that it would ring when the result is in unwanted range***
- 4.When the measurement fails due to some factors, it can show relative error info.***
- 5.Integrated with SPO2 function, SPO2 probe is optional***
- 6.Able to communicate with PC, you can make review, analyses, trend graph and report printing with the software.***

Considerateness

- 1.You can set time so that you can know the time for every measurement, and the clocking keeps going when you change the battery in a short time***
- 2.You can shift between the measurement units of mmHg and kPa according to your own habit.***
- 3.It can remind you when power is low, and can be power off automatically when left idle up to 5 minutes.***

Technical Specs:

Measurement

Measurement Method: Oscillometry

Measurement Position: Upper Arm

Measurement Range:

Pressure:0~290mmHg(0~36.0kPa),Pulse

Rate:40~240bpm,Spo2:35%~100%

Measurement Accuracy: Blood Pressure: ± 3 mmHg, Pulse

Rate: large of 3bpm

and $\pm 3\%$,Spo2: $\pm 2\%$ in 70%~100% range

Display

2.8" colour LCD

Power Supply

4 AA Size Alkaline Battery or optional AC/DC Adapter

Configuration

Standard Conf.: Adult cuff, software disc, user manual, USB

line, 4 AA Size

alkaline battery

Features Briefing:

1.One touch auto inflate

2.Colorful, high resolution, large LCD

3.Three-user data system,99 records for each user


- 4. Mode for adult, pediatric and neonatal selectable*
- 5. Alarm setting*
- 6. Error detection and showing*
- 7. Optional SPO2 function*
- 8. Communication with PC*
- 9. Time set, Clocking undisturbed in short time battery changing*
- 10. mmHg and kPa selectable*
- 11. power indication, auto power off when left idle*

Software Functions:


- 1, Upload data to computer*
- 2, Show trend graph, analytical diagrams*
- 3, Edit add notes to every piece of data*
- 4, Edit examinees' information, medical advices, condition descriptions, present drug usage info etc.*
- 5, Print report, support print preview*

Standard Accessories :

- Main units : 1 pc*
- Adult Cuff: 1 pc*
- Software CD: 1pc*
- USB cable :1 pc*
- User manual: 1 pc*

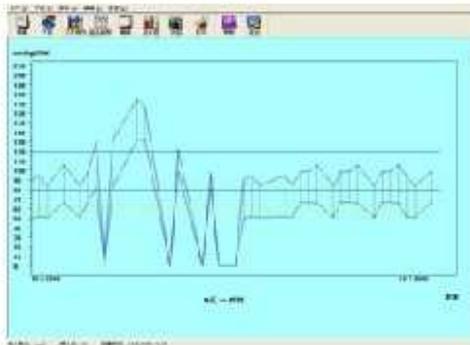


Trend

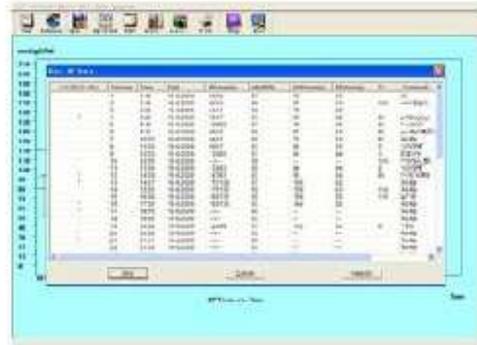

Editdata



Pie

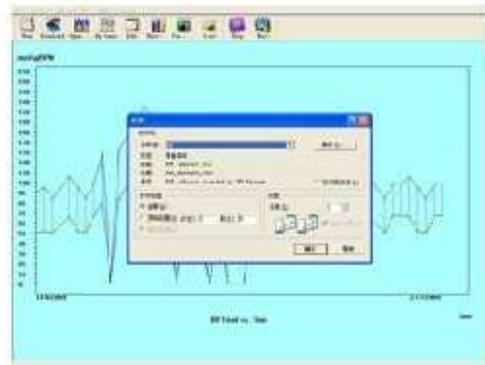



Print report

Multi-interface





Bar Graph for Pressure and Pulse Rate Analysis

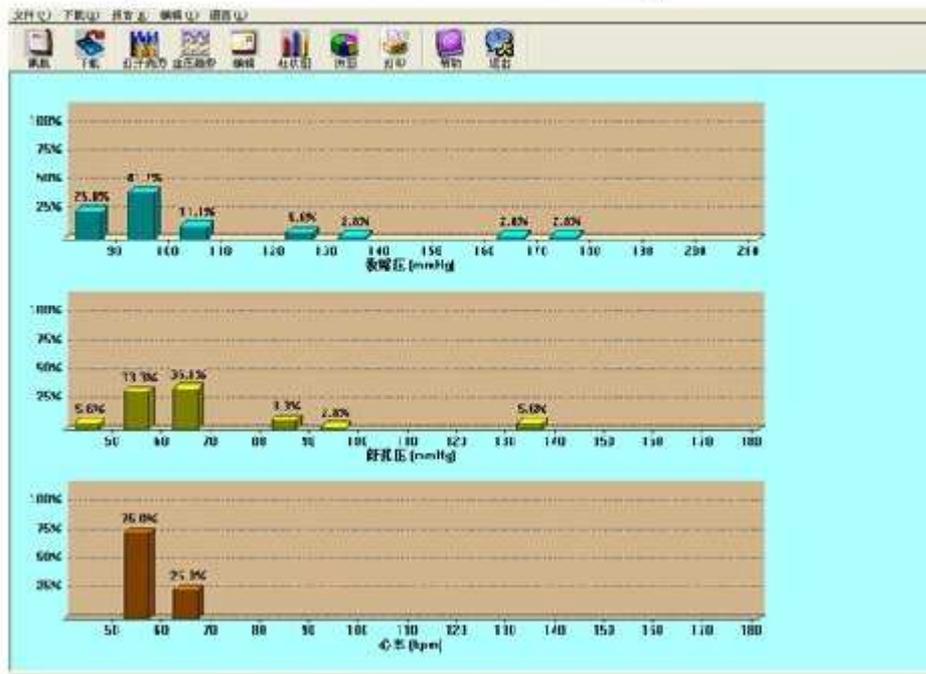

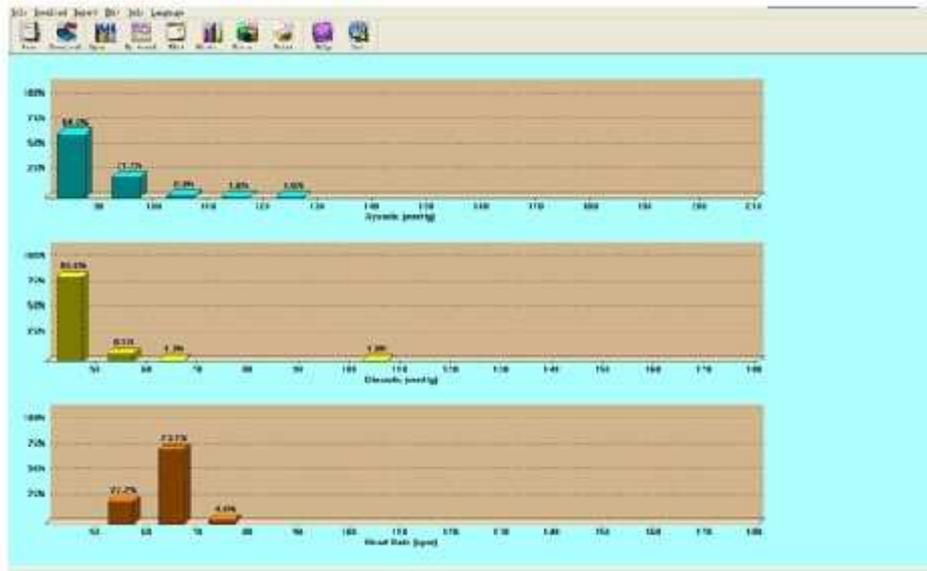

Pie Graph for Pressure and Pulse Rate Analysis

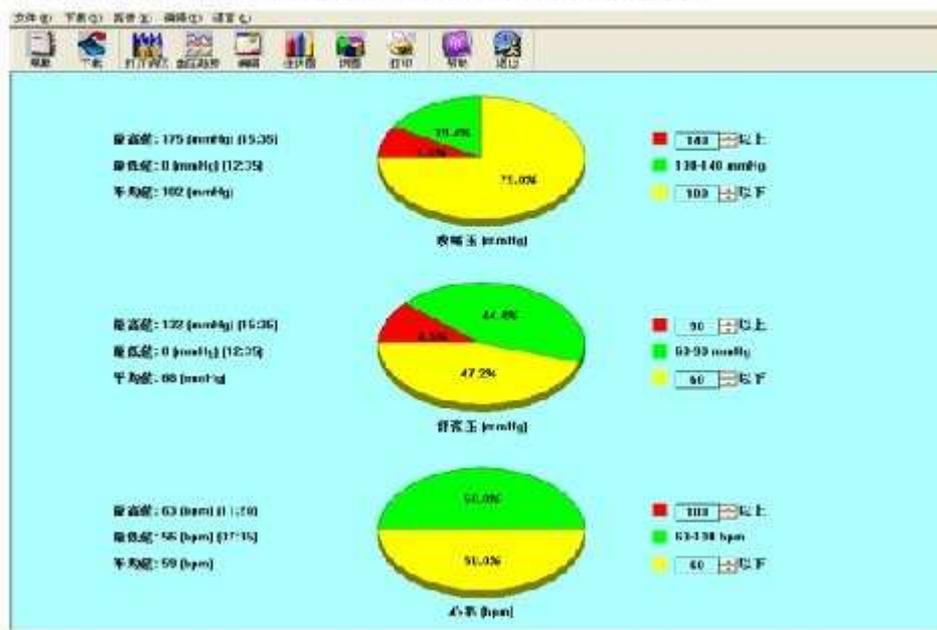

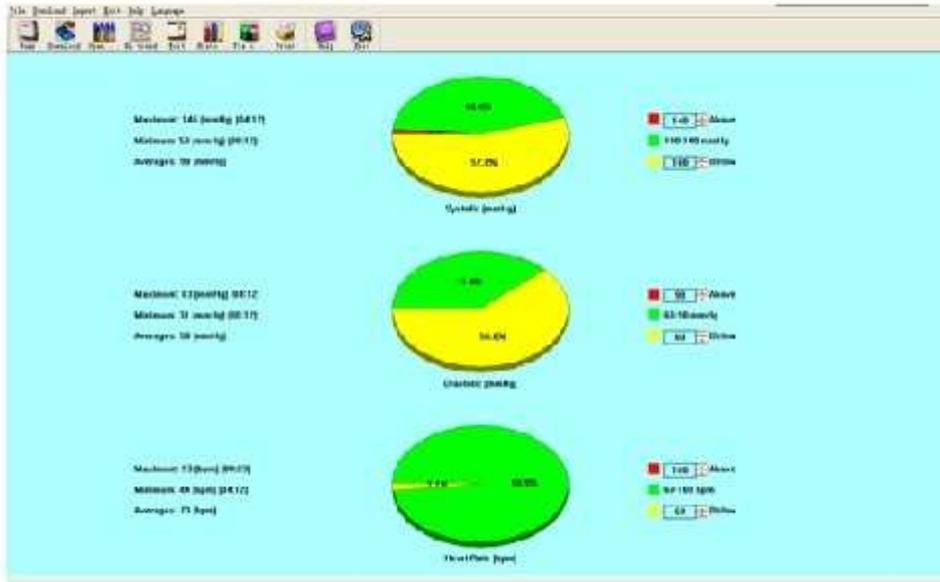

Pressure Trend Graph

