


# BLOOD AND INFUSION WARMER

## FT70

### Economical fluid warmer for:

- Routine infusion in slow or medium rate
- Enteral and parenteral nutrition and other solutions
- Infusion to children and neonates
- Normal infusion in cold conditions

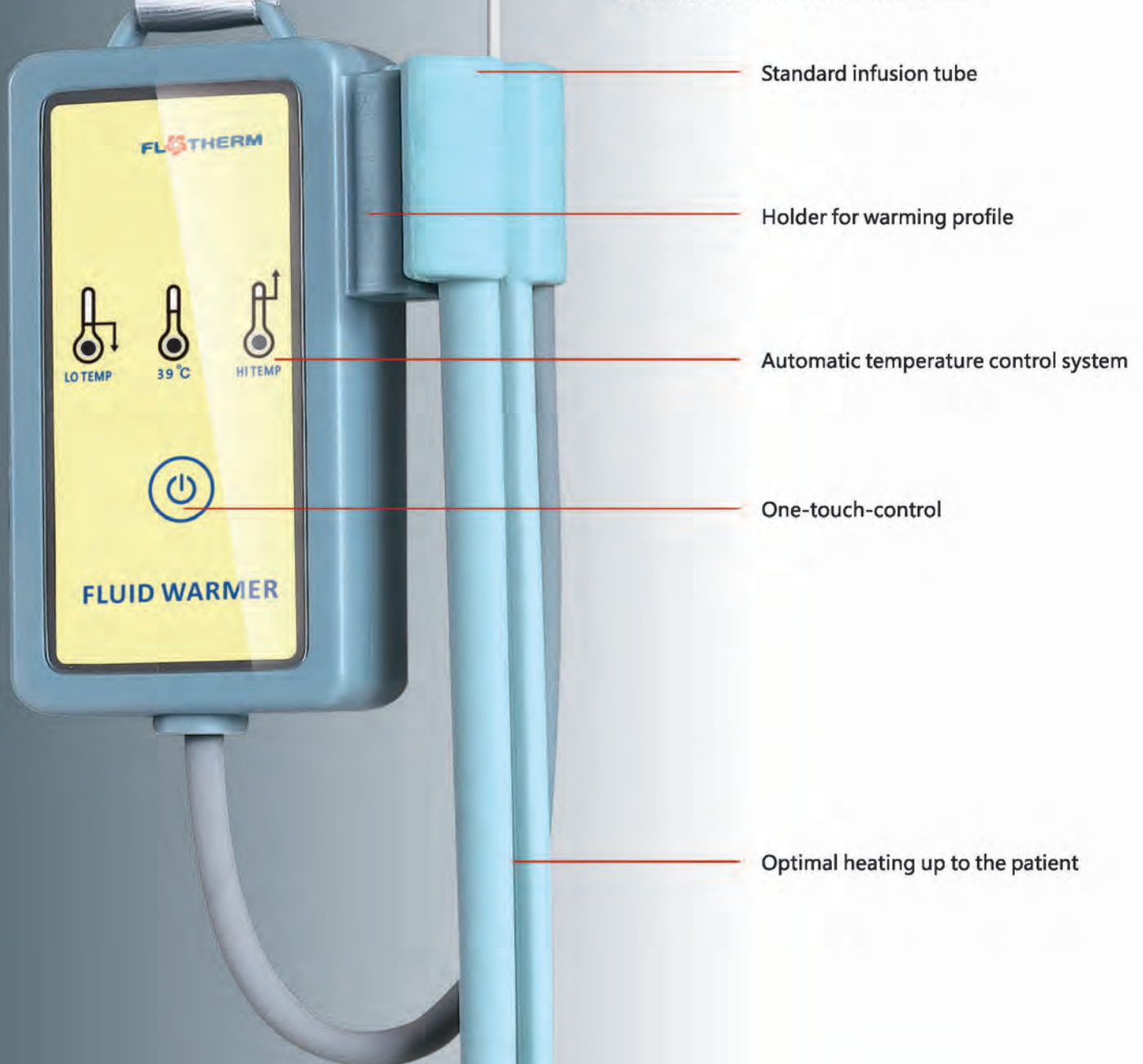

# BLOOD AND INFUSION WARMER

## FT70

### FEATURES

#### Open system:

- To use independently or be equipped to other warming system, infusion pumps or feeding pump
- Accepts standard IV tube, no special disposables needed

#### Heating up to the patient:

- IV tube is completely wrapped in, no heat loss

#### One-touch-control:

- Temperature automatically controlled
- Easy to operate
- Effective and gentle heating
- Compact in size

### TECHNICAL SPECIFICATIONS

Model:	FT70
Temperature setting:	39°C ex-works
Power supply:	a.c.100~240V/50~60Hz
Power consumption:	Max. 50VA
Type of protection against electric shock:	Class I
Degree of protection against electric shock:	BF Applied part; Defibrillation-protected
Degree of protection against ingress of liquids:	IPX2
Temperature accuracy:	±2.0°C
Overheat protection:	42°C
Low temperature alarm:	37°C
Warming up time:	From 20°C to 39°C approx.3 min.
Operating mode:	Continuous
Dimension (W*D*H):	70×25×90mm
Net Weight:	0.7kg
Warming profile:	length 1000mm; compatible to IV tube 3.5-5.0mm O.D.

\* Subject to technical modifications without notice

**keewell**  
Medical Technology

Keewell Medical Technology Co.,Ltd.

Xiananyi Industrial Zone, Pingzhou, Nanhai, Foshan,  
Guangdong 528251, P.R.China

T.: +86 757 86260839 F.: +86 757 86260100  
marketing@keewell.cn http://en.keewell.cn

CE 0482

Distributor: