

NEW HOSPIVAC

The NEW HOSPIVAC Line has been designed for operating theatres for the aspiration of body liquids and for applications in gynecology and dermatology (liposuction). Electronically managed footswitch control and jar selector available upon request.

Convenient Handle
for Easy Movement

Compartment for
Accessories

Foot-Switch Operation
Electronic Button,
Double Mode ("FS" and
"FULL" Models Only)

Vacuum Indicator
(kPa and Bar)

Antibacterial and
Hydrophobic Filter

Jar Selector One-Touch
("FULL" Models Only)

Vacuum Regulator


Duble Suction Inlet
("Full" model Only)

Collection Jar in
Autoclavable PC with
Overflow Valve System

Silicone Autoclavable
Tubes

Four Antistatic Castors
(Two with Brake
System)

CE 0123

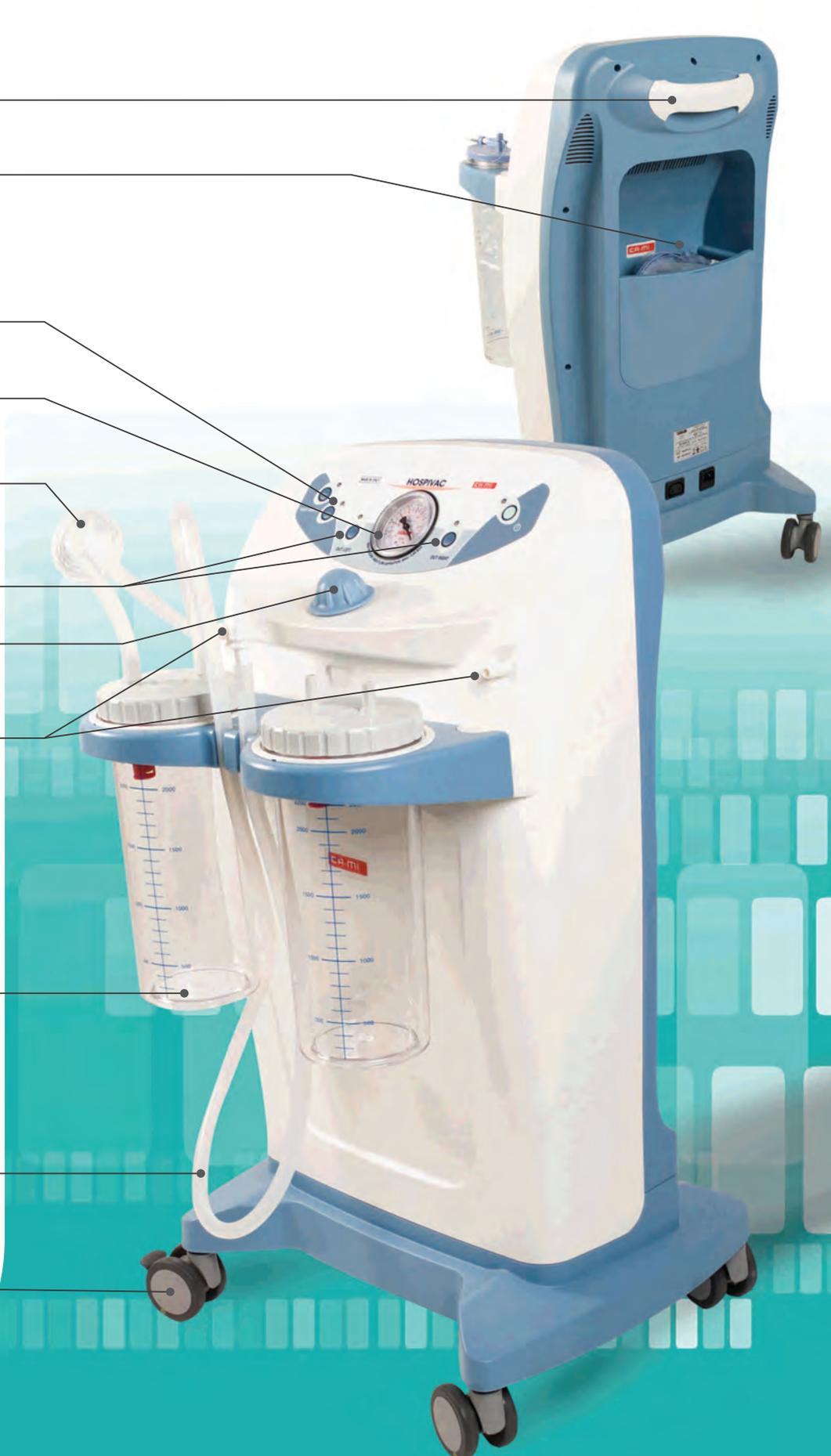

NEW HOSPIVAC 350



TECHNICAL FEATURES

Motor	Oiless and maintenance-free piston pump
Power Feeding	230V-50Hz Other configurations available upon request
ISO 10079-1 Classification	HIGH VACUUM / HIGH FLOW
Max Vacuum (adjustable)	-0,90 Bar -90 kPa -675 mmHg
Max free air flow rate	60 l/min
Noise Level	51,7 dB
Power consumption	230 VA
Fuse	1 x F4 A 250V
Duty cycle	Non-stop operation
Weight	Kg 13
Size	cm 46 x 85 x 42
Years of Warranty	2
Shipping carton	1
Place of Manufacturing	Italy

MAIN APPLICATIONS

Emergency Dept.

Operating Theatre

General Surgery

Neurosurgery

Liposuction

Dental Practice

Gynecology

NEW HOSPIVAC 400



TECHNICAL FEATURES

Motor	Oiless and maintenance-free piston pump
Power Feeding	230V-50Hz Other configurations available upon request
ISO 10079-1 Classification	HIGH VACUUM / HIGH FLOW
Max Vacuum (adjustable)	-0,90 Bar -90 kPa -675 mmHg
Max free air flow rate	90 l/min
Noise Level	46,4 dB
Power consumption	300 VA
Fuse	1 x F4 A 250V
Duty cycle	Non-stop operation
Weight	Kg 20
Size	cm 46 x 85 x 42
Years of Warranty	2
Shipping carton	1
Place of Manufacturing	Italy

MAIN APPLICATIONS

Emergency Dept.

Operating Theatre

General Surgery

Neurosurgery

Liposuction

Intensive Care Unit

Gynecology

11