

MIMSAL



Healthcare Lighting

2021

MIMSAL

MIMSAL

CATALOGUE HEALTHCARE LIGHTING 2021

© MIMSAL TRADE S.L.

C. Mollet 17, Pol. Ind. Palou Nord
08401 Granollers (Barcelona) - Spain
Tel. +34 930 139 860
mimsal@mimsal.com
www.mimsal.com



MIMSAL

MIMSAL is a specialized lighting solutions supplier for the healthcare sector. We offer a wide range of luminaires designed to satisfy the diversity of the end-user needs and for each activity sector, by offering technical advice and commercial support to the professionals / prescribers of the healthcare.

Our products are engineered for easy mounting, with modern electronic components and light sources for the best energy efficiency and economy.

Our products and solutions are manufactured and certified in accordance with Regulation (EU) 2017/745 of the European Parliament and of the Council on medical devices.

MIMSAL complies with the quality certification ISO 9001 and ISO 13485.



LUPA LED H.F

LED magnifier with round lens, Lupa LED H.F is intended when performing any meticulous work in need of great accuracy. It is ideal for hospitals, measured studies, laboratories, aesthetic clinics, etc.

Functionality and Performance

Generates a light without dark spots.
Internal spring balanced arms that allow different positions.
Switch on the head.
Round head with convex lens, adjustable and simple mechanism

Mounting Options

The standard version includes AH table clamp.

As for accessories, we offer a wide selection of supports such as wall supports, wall extension arms, trolleys stand and rail brackets.

Cleaning and Safety

Head exempt of heat.
Protective cover on the head with hinge mechanism to prevent from dust and scratches that may damage the lens.

Economy

High performance, low power 6W.
Long life, 40,000 hours.
Durable materials.
Very low maintenance costs.



LUPA LED H.F

Technical Details		(Cod. 09340)
Light source	LED H.F 6W	
Lamp lifetime	40.000 hours	
Colour temperature	6.500 °K	
Glass lens	5 diopters	
Lens diameter	127 mm	
Arm length	90 cm	
Weight	3 kg	
Colour: White	White	



Trolley stand
8,8 kg Cod. 09651





Wall support
Cod. B0101800




AH table clamp
Cod. B010190



Rail Plus bracket
Cod. 82026




Extension wall arm
Cod. B010150

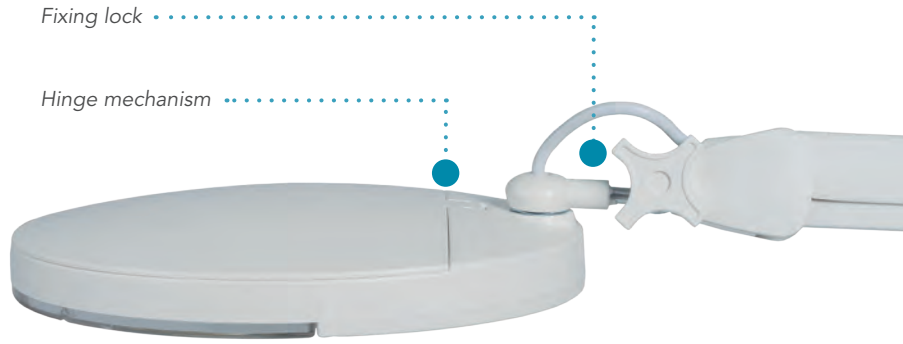


Fixing lock

Hinge mechanism

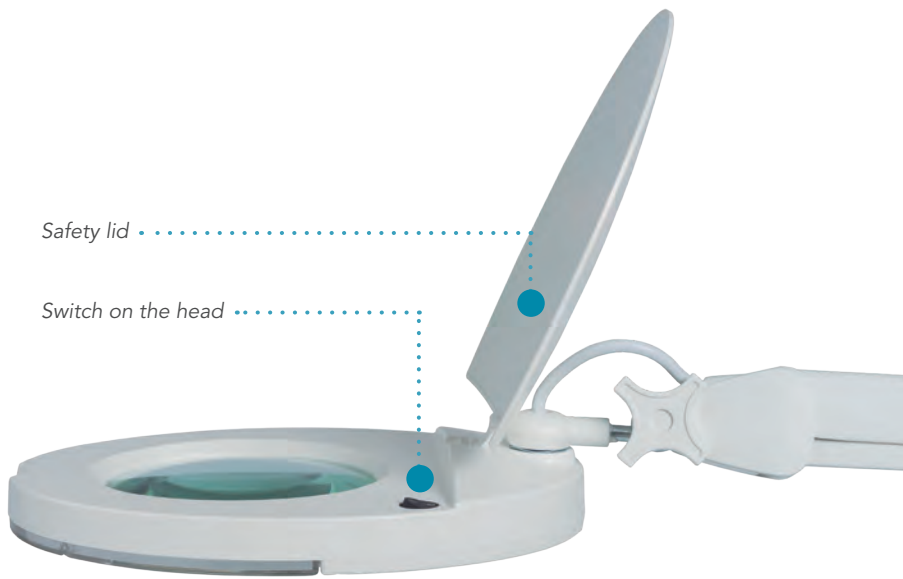

Safety lid

Switch on the head

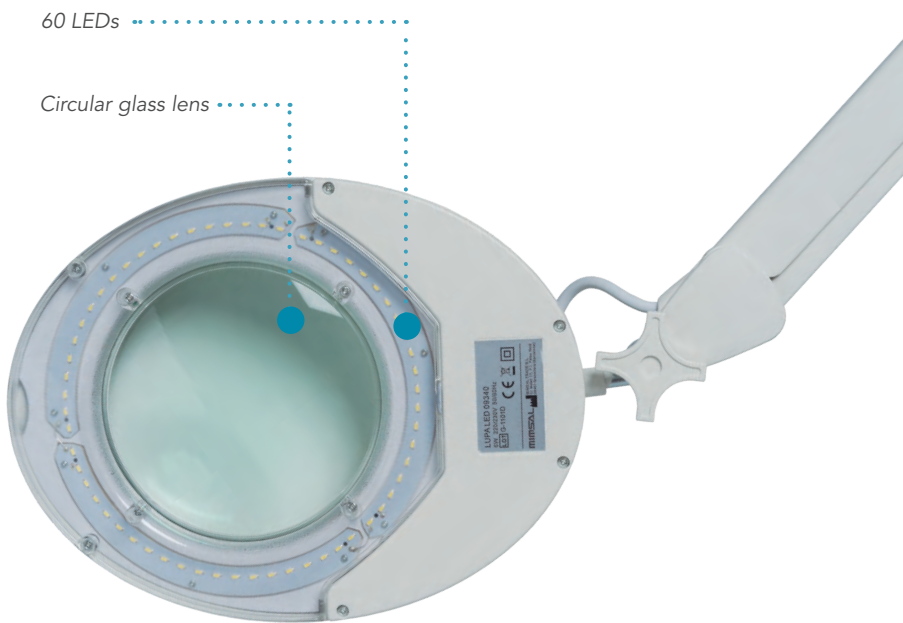

60 LEDs

Circular glass lens

LUPA H.F

Cold light luminaire with magnifying lens, essential for high precision work. An ideal use for hospitals, medical centers, laboratories and aesthetic clinics.

Light Source

Low consumption 22W circular fluorescent tube.

Produces a light without dark spots and does not overheat.

Lifetime of 10.000 hours.

Functionality and performances

Adjustable height articulated arm that allows easy movement in different positions.

Head with biconvex circular lens of 5 diopters, adjustable and simple mechanism.

Switch on the head of the luminaire.

Security and hygiene

Head with protective cover to prevent dust or scratches that may damage the lens.

Mounting Options

Standard version includes AH table clamp.

As for accessories, we offer a wide selection of supports such as wall supports, wall extension arms, trolleys stand and rail brackets.



LUPA H.F

Technical Details	(Cod. 09332)
Light source	FL 22W
Lamp lifetime	10.000 hours
Colour temperature	6.500 °K
Glass lens	5 diopters
Lens diameter	125 mm
Arm length	90 cm
Weight	3 kg
Colour	White

Trolley stand
8,8 kg Cod. 09651





Wall support
Cod. B0101800



AH table clamp
Cod. B010190



Rail Plus bracket
Cod. 82026



Extension wall arm
Cod. B0101S0



ACCESSORIES

WALL B SUPPORT

Wall support.
Separation between axis-wal: 7cm
Cod. B0101800




AH TABLE CLAMP

Table fixing.
Cod. B010190



RAIL PLUS BRACKET

To adapt to U.C.I rail.
Cod. 82026

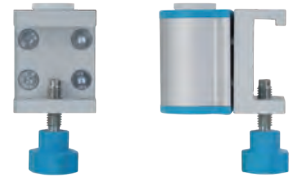

TROLLEY STAND

H=92 cm Ø 61 cm.
4 kg Cod. 09655
8,8 kg Cod. 09651



EXTENSION ARM + WALL B SUPPORT

It provides 40 cm rotation from the wall.
Cod. B010150





MIMSAL

MIMSAL TRADE S.L.

C. Mollet 17, Pol. Ind. Palou Nord
08401 Granollers (Barcelona) - Spain
mimsal@mimsal.com - www.mimsal.com