



Save • Someone's **life**



Defi[®] 5S +

Automated External Defibrillator
Fully Automatic

Defi[®] 5S +

Automated External Defibrillator

ONE BUTTON FOR **HEART**..



Start the device

1

Turn on the device.

The voice prompt immediately begins with instructions on how to operate the AED.



Attach the pads to the patient

2

After the electrodes are attached,

the AED immediately begins analysing the ECG.

The device does not proceed to the next step until you have properly attached the electrodes.



Select energy

3

International resuscitation guidelines recommend reducing energy for children to avoid burning the heart muscle. The Meditech Defi 5S series can reduce the energy to 50.200J by simply selecting the energy range.



Automatic Shock

4

When the shock is needed,


The device will deliver the electric shock to the patient.

WHY

Advantages




1 LED Indicators
Color LED indicators for each step

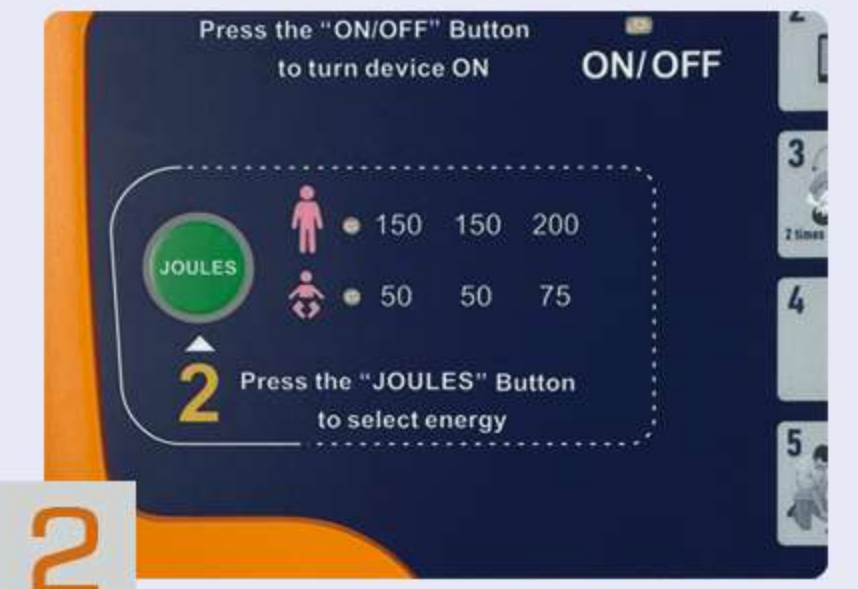

2 Energy Sequence Selectable
Energy sequence selectable 50J to 200J for different age groups and situations




3 Loud Speaker
High quality speaker design Provides great and loud sound that makes it easy to hear the device's instructions in public places

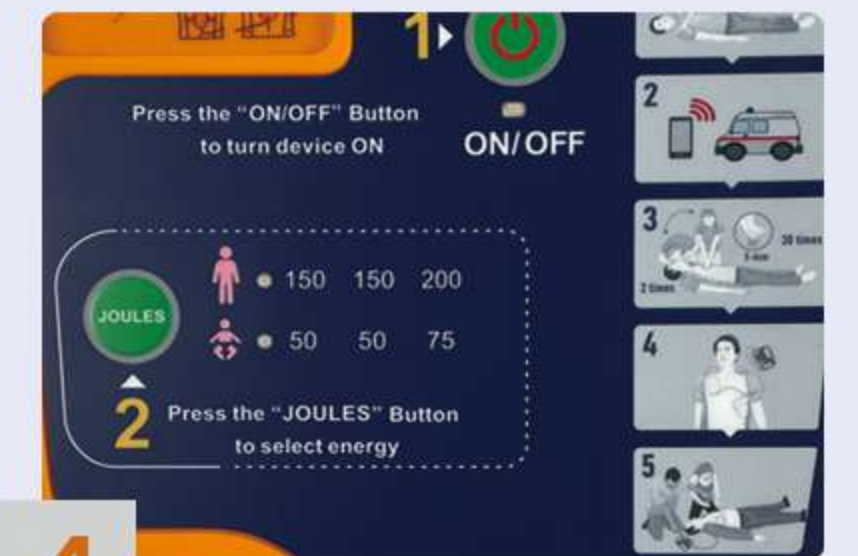

4 Guide line
Visual steps for using an Automated External Defibrillator



5 Written use Instruction
Written instructions on the device for each step and key

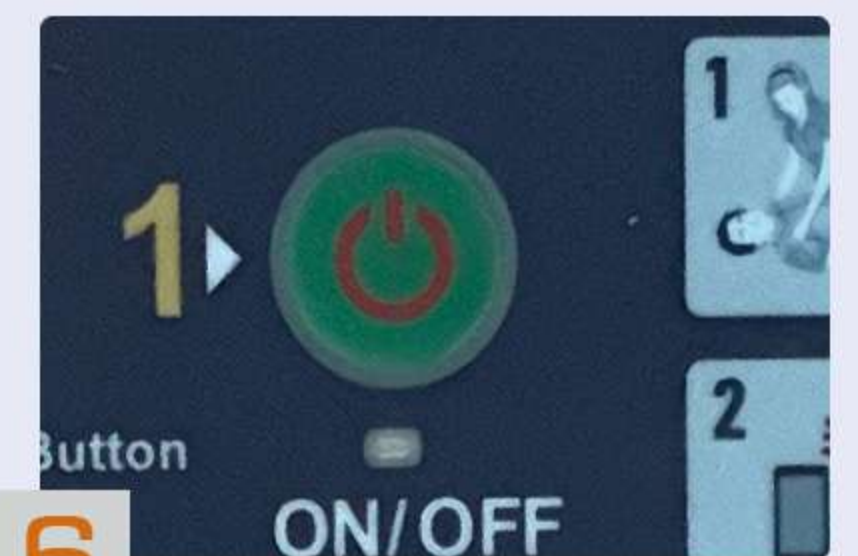

6 Buttons
Waterproof buttons and soft



2 Languages
Selectable



Child Mode
50.50.75J



Adult Mode
150.150.200J



Infra-Red
Data transfer



Self test
Day.Week.Month



IP55
Water and dust resistance



Sound volume
5 Levels



Carry Case
for AED and first Aid



Built-in Memory
ECG waveform



Biphasic function

- Fully automatic system
- 5 volume levels (by the end user)
- Button for adult/child mode
- Use of AED in pocket (does not need to be opened)
- Extra pocket in the bag for first aid products
- IPX55 water and dust resistant
- 2 languages to choose from (by end user)
- Infrared for data transmission (standard)
- Daily, weekly and monthly self test

Data Communications

The AED memory software assists the AED in receiving data sent via infrared. The Meditech AED Memory Software supports reading and setting ECG waves, setting the AED language, reading, analysing, and deleting memory data can also be performed with the Meditech AED Memory Software.

Case Management

Meditech's AED storage software makes the selected data file a case record and stores cases in the database ACCESS. Case management supports importing and deleting data that can respond to the ECG waveform of the data file by selecting cases.

Memory Playback



The interface displays status information recorded during AED operation, including type, event time, description, and waves. Event types include: Power On Event, Self Test Event, Analysis Event, Charge Event, Discharge Event, CPR start event, CPR stop event. The event time displayed indicates when the event occurred. The description displays the parameter information. "OHM" represents the detected impedance of the patient. The parameter "Ppower" stands for the voltage of the completed charge; "Shock required" means that a shock is recommended after the analysis is complete "No shock required" means that no shock is recommended after the analysis.

ECG

Waveform gain can be set to 5mm/mv, 10mm/mv, 20mm/mv; Waveform speed can be set to 12.5 mm/s, 25 mm/s, 50 mm/s;


Print

Can preview the selected page and print basic information

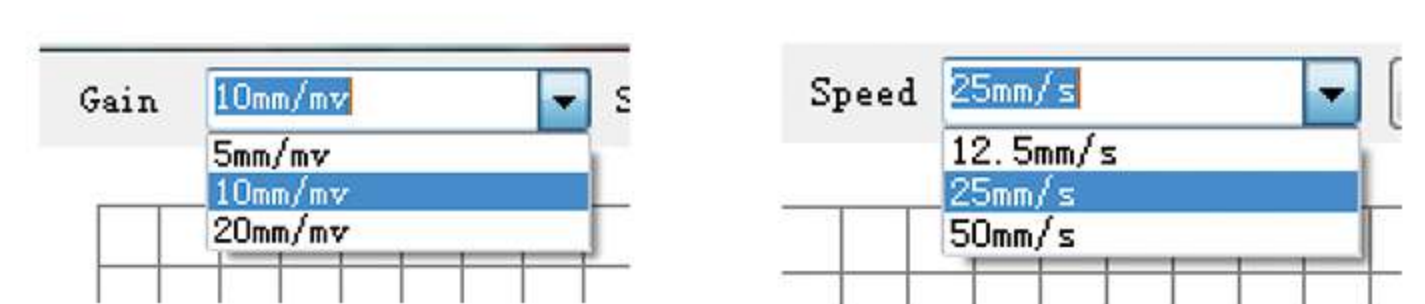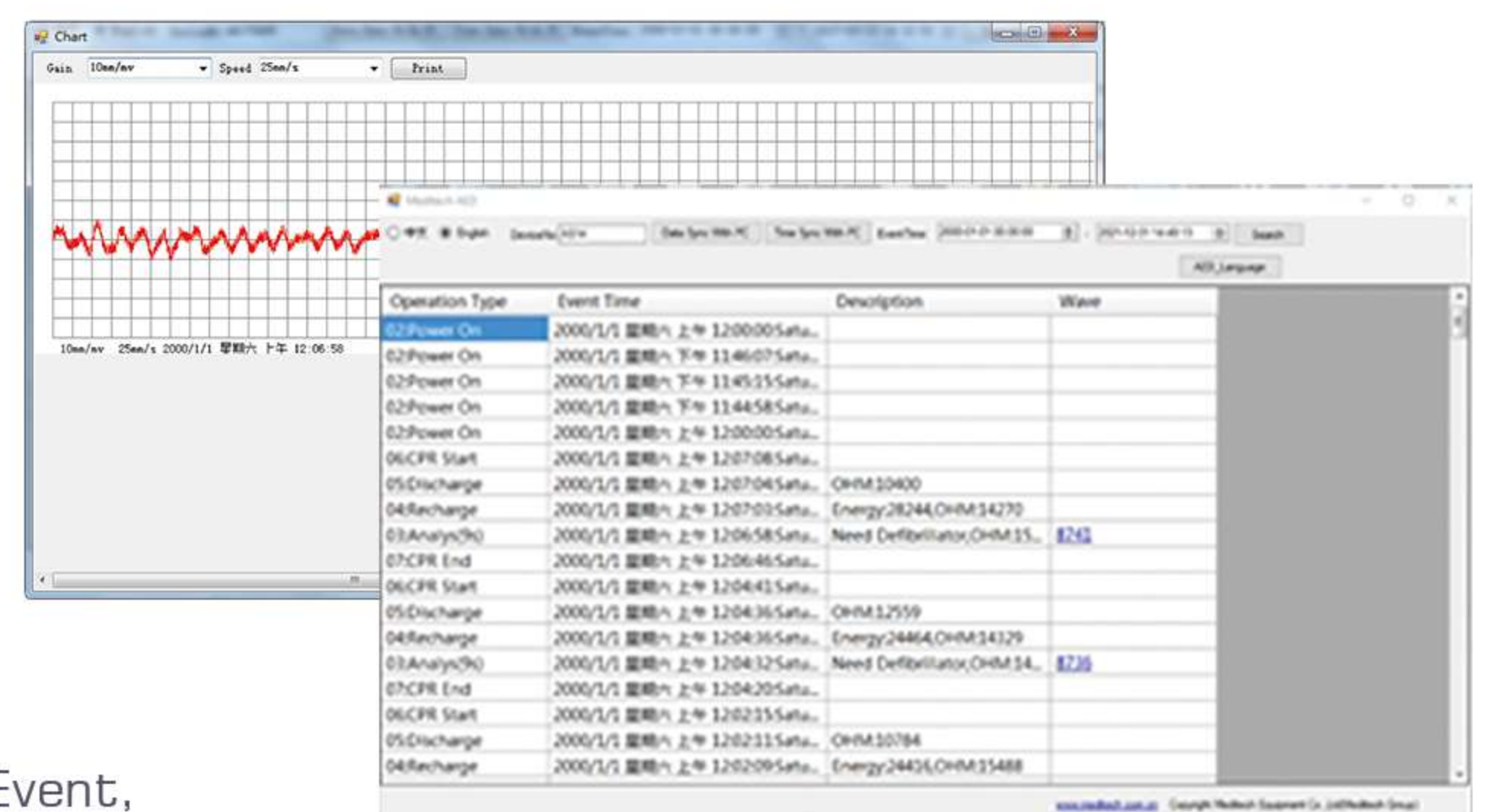

Technical Specifications

Self-test Mode	Daily, Weekly, Monthly
Waveform	Adult, Child
Energy	Biphasic truncated exponential
Energy Sequence	200-Joules max.
	Child Mode: 50.50.75 Joules
	Adult Mode: 150.150.200 Joules
Charge time (New, at 25°C)	Less than 8 sec. to 150J Less than 12 sec. to 200J
Voice Prompt	Extensive voice prompt
Visual Indicators	LED prompts
Control	Two buttons: ON/OFF, Mode
ECG storage	1,800 events
Data transmission	Infrared
Power Type	Battery non-rechargeable 12V, 2800mAh Li-MnO2 Cell
Operating Temperature	0°C to 40°C
Operating Humidity	Relative humidity between 30% and 95% (non-condensing)
Storage Temperature (without battery)	-20°C to 55°C
Storage Humidity (without battery)	Up to 93% (non-condensing) Physical
Storage Temperature (Disposable Pad)	0~50C
Dimensions	Depth (D): 80mm Width (W): 240mm Length (L): 300mm
Weight (with battery)	1.9Kg (4.2 Pounds)

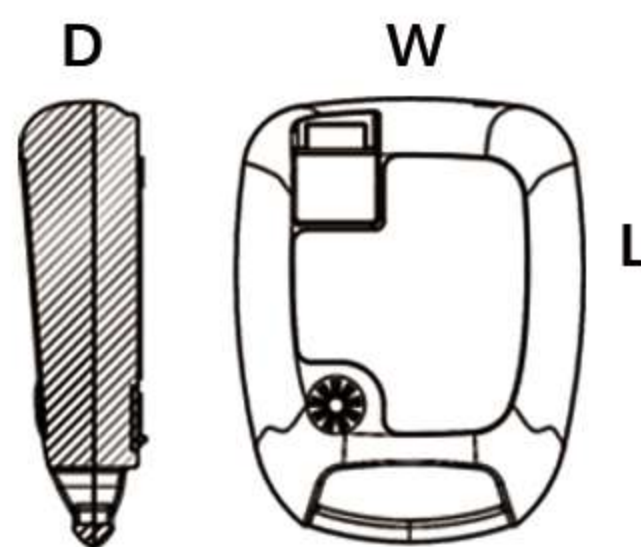

Accessories

Battery

CR123A - 4x2 2.8Ah 12Volt non-rechargeable battery

Disposable Pad



Designed for external defibrillation, non-invasive pacing, synchronized cardioversion, and monitoring of adults and children weighing more than 20 kg (55 lbs).

Carry Case



- 1 Reflective sign
- 2 Soft Hand
- 3 Speaker protection
- 4 Transparent panel
- 5 Socket Protector





MEDITECH[®]
EQ

Meditech Equipment Co.,Ltd (Meditech Group)
Address: 140 Zhuzhou Rd, Qingdao ,
Shandong Province, P. R. C
Tel: (+86) 532 85832673
Fax: (+86) 532 81705331
Email: sales@meditech.cn
www.meditech.com.cn

Europe: (+44) 020 3239 3734
Email: UK@meditech.cn

USA : (+1) 509 588 7272
Email: US@meditech.cn

Middle East: (+2) 02 27740444
Email: Eg@meditech.cn

CE 0482

