

MEDITECH[®]
ED

Advanced Medical Technology




Advantage

Aあ

2 Languages
selectable



Child Mode

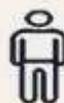

Adult Mode



Program
4 training modes



Sound volume
5 Levels



Carry Case
for AED and first Aid

Defi[®] 5S Trainer AED Trainer



FOUR
Different scenarios



Features

Defi[®] 5S Trainer

- 15 languages to choose from (by end user)
- 3 volume levels (by the end user)
- Button for adult/child mode
- 4 training modes
- Use of the trainer in the pocket (does not need to be opened)
- Extra compartment in the bag for first aid products
- Extensive voice prompts for the operator
- Manually switchable language
- Runs on standard AA x 4 batteries

Why are AED training important?



AED training is important because it equips individuals with the knowledge and skills necessary to save a life in the event of sudden cardiac arrest. With proper training, individuals can recognize the signs of cardiac arrest and use an AED to provide life-saving support until emergency services arrive.



How to use



- 1 Start the device**
 Turn on the power.
 The voice instruction starts immediately for guide on how to operate AED.

- 2 Attach the pads**
 After you connect the pads, the trainer starts working with the selected mode (one of four).



- 3 Selection mode**

Mode 1: Shock indication and CPR voice instruction

Mode 2: No shock indication and CPR voice instruction

Mode 3: Shock indication and no CPR voice instruction

Mode 4: No shock indication and no CPR voice instruction

- 4 Press the button**
 When the trainer asks you to give a shock (the trainer does not have an electric shock). Press the red shock button





- 1 Reflective sign
- 2 Soft Hand
- 3 Speaker protection
- 4 Transparent panel
- 5 Socket Protector



Specifications ■ ■ ■


Mode: Adult, Child

Dimensions: Depth (D): 80mm

Width (W): 240mm

Length (L): 300mm

Weight (with battery): 1.9Kg (4.2 Pounds)

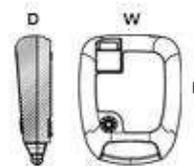

■ Battery

4 x AA Battery (6Volt)

■ Reusable pads



Language

Defi[®] 5S Trainer

More than 45 languages are available



Speaks
your
language



Fast action saves lives!

