

HEINE BETA X Otoscope

The world's first otoscope with an additional integrated examination light and exchangeable optical modules.

Large diameter viewing window

for extra-large viewing area

The new HEINE X CHANGE System

Expand the function of your otoscope with easily interchangeable X CHANGE modules:

► HEINE X CHANGE System 3x Magnifier

- 3x magnification
- Perfect all-rounder magnification option
- Multi-coated precision optics for crisp images
- Quick and easy view with a large field of view

► HEINE X CHANGE System 4.2x Magnifier

- 4.2x magnification
- Advanced multi-lens optical system for higher magnification and resolution
- Excellent depth of view with advanced spatial recognition to easily identify even the smallest anatomical structures and foreign bodies
- Crystal clear images no matter if you prefer a closer (traditional) or greater working distance from the instrument

► HEINE X CHANGE System Smartphone Adapter

- Easily capture digital images with any smartphone

More HEINE X CHANGE System modules coming soon!

Register at heine.com/xchange to be the first to receive updates

Integrated insufflation port

for tympanic mobility testing

Intuitive illumination control functions

- Easily switch between otoscope and the exclusive HEINE inSPECT illumination options
- Select from two light intensities for both illumination options
- Simple one-finger operation

Smart power source options

- Charge status indicator
- Automatic switch-off function
- Rechargeable with USB-C or charging station



HEINE BETA X Otoscope

HEINE LED^{HQ} and fiber optic technology for optimal illumination

- Homogeneous, very bright illumination with true colour rendering for high contrast and extremely detailed views
- Fiber optic provides an unobstructed view of the ear canal and tympanum
- Extended battery operating time and a virtually unlimited operating life of the LED^{HQ}

Compatible with EcoTips

HEINE AllSpec disposable tips made from recycled plastic

smooth SURFACE

Safe, hygienic reprocessing thanks to smoothSURFACE and recessed connections




Product data

The new HEINE inSPECT additional examination light for the eyes, throat and more

- A patent pending world exclusive: medically tested and safe for examining the anterior ocular media and the pupillary light reflex
- Ideal for illuminating deep anatomical structures in the mouth and throat

Durable, lightweight and ergonomic design

- Solid, integrated aluminium frame
- Housing is thermo-moulded out of a special high-impact material blend optimised for hygienic reprocessing
- Precisely balanced for excellent ergonomics and secure handling

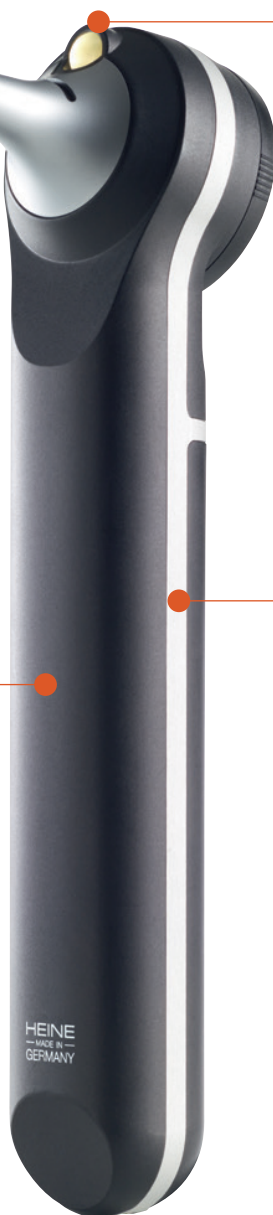

HEINE Optotechnik GmbH & Co. KG

Dornierstr. 6, 82205 Gilching, Germany
www.heine.com



We reserve the right to change specification without notice.

04|25. A-000.00.273 EN

QUALITY
MADE IN GERMANY

