

# KANMED<sup>o</sup> WARMCLOUD


## PATIENT WARMING SYSTEM



Efficient patient warming system with pressure sore prevention




Optimal warming and pressure relief

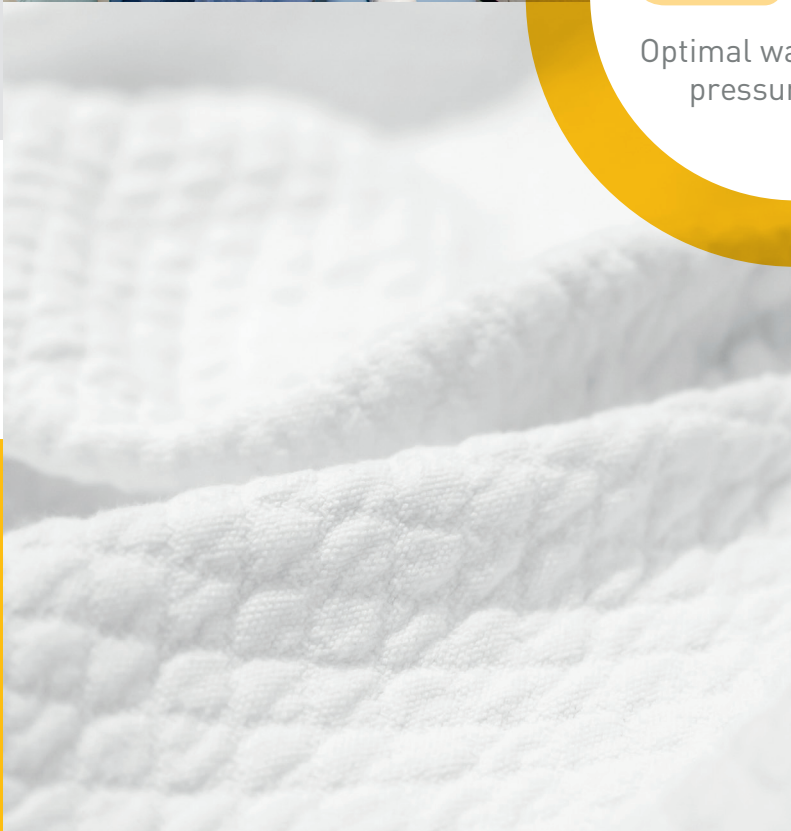

#### KANMED WARMCLOUD

ensures efficient and thorough warming of the patient during long surgical procedures. The soft, disposable, single patient use mattress, reduces the risk of pressure sores – especially during very long operations.

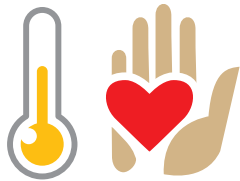

## KANMED<sup>o</sup> WARMCLOUD

# Patient Warming System with Pressure Sore Prevention

Kanmed WarmCloud is a unique system that ensures patients stay warm during surgery and reduces the risk of pressure sores. The Kanmed WarmCloud system consists of a single-use warm air mattress, a main unit and an optional remote control for operating the system. The warm air mattress placed under the patient, lifts the patient creating a soft, comfortable and temperature-controlled environment.

#### **The best solution for the patient**

The system makes the cold, hard operating table warm and comfortable in no time at all. This helps the patient relax. It dramatically reduces the risk of pressure sores, especially during very long procedures.

#### **The best solution for the anesthesia staff**

Place the patient on the mattress and start the unit. Very quickly, your patient will be warm, comfortable and cooperative. WarmCloud can follow the patient from preparation room to operation room and into post-op/ICU

#### **The best solution for the surgeon**

Because the mattress is under the patient the surgeons will have easy access to surgical sites. The design prevents warming of the operating staff. It also prevents disturbance of the laminar air flow system. The WarmCloud is very quiet.


#### **The best solution for the hospital**

Efficient patient warming shortens postoperative stay, since the well-known hypothermia-induced complications are eliminated. Reduction of complicated procedures for pressure sore prevention saves time and costs. The result is satisfied patients.

Superior patient comfort

Closed circuit. No air contamination.



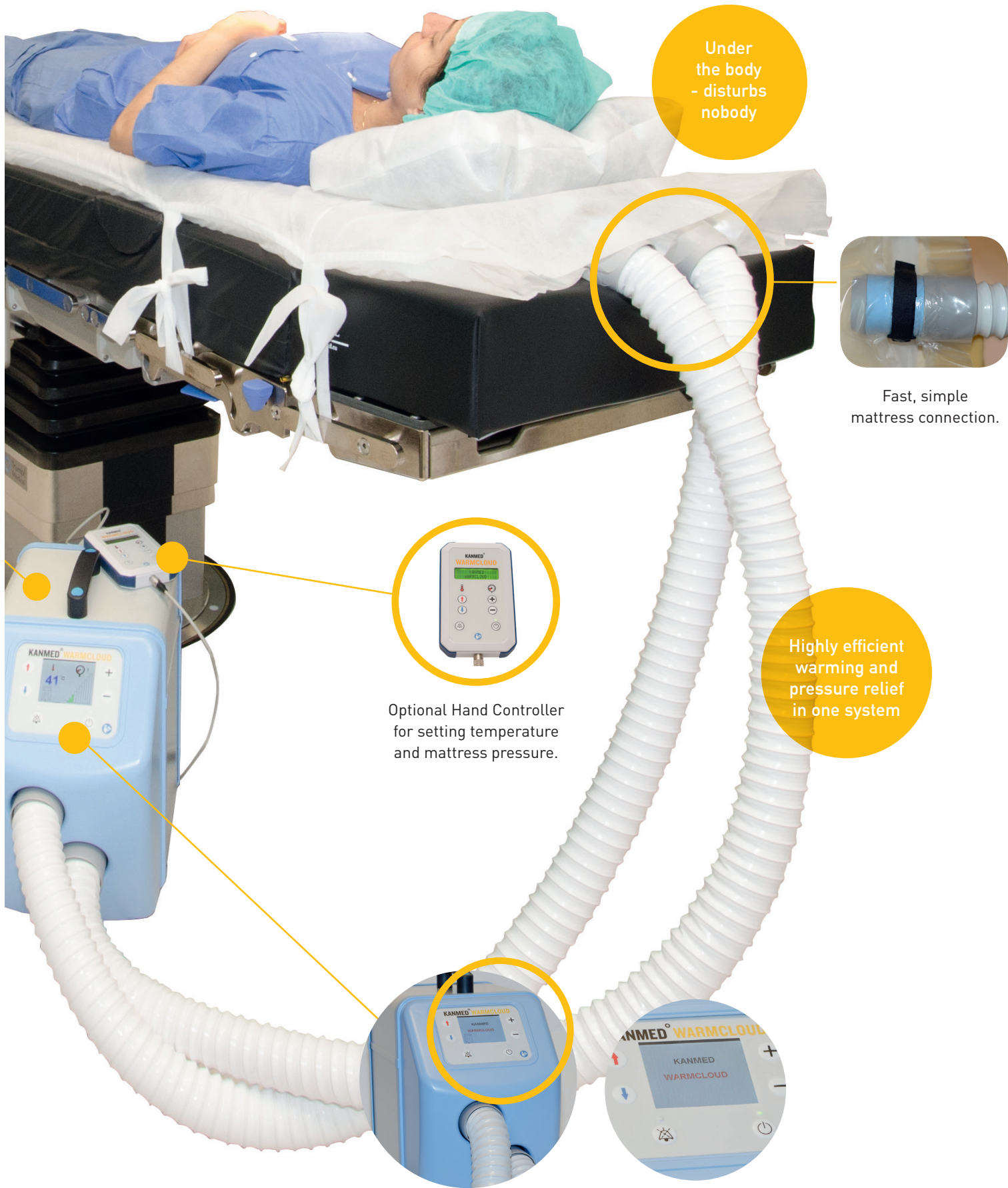

Under the body - disturbs nobody



Fast, simple mattress connection.



Optional Hand Controller for setting temperature and mattress pressure.

Highly efficient warming and pressure relief in one system



The quiet Control Unit can be placed under the head or foot end of the operating table. It may also be hung elsewhere using its handle and an adapter. A wheel base with tube support is available.

# KANMED<sup>o</sup> WARMCLOUD

## COMPLETE SYSTEM WITH STANDARD ACCESSORIES

Art.no 9000 230VAC, 50/60Hz. (115VAC, 50/60Hz is also available).

### Main Unit

Dimensions: Width 295 mm. Height 370 mm. Length 210 mm

Weight: About 12 kg incl. wheelbase

Length of Air Hoses: 150 cm

Temperature range: 31°C–42°C

Pressure range: 1–10 (about 5–50 mb)

Length of cable: 2,9 m

### Remote control (option)

Dimensions / weight: Width 90 mm. Height 25 mm. Length 150 mm. Weight 200g

## WARM-AIR MATTRESS, FULL LENGTH

Art no. OT-600-222

Dimensions: Length 200 cm. Width 55 cm

Weight: 0,95 kg

Package 10 pc. / box

Material: EVA plastic with non-woven cover, 5 side straps

## WARM-AIR MATTRESS, SHORT

Art no. OT-600-211

Dimensions: Length 110 cm. Width 50 cm

Weight: 0,45 kg

Package 10 pc. / box

Material: EVA plastic with non-woven cover, 4 side straps

- The mattress is a single-use item. Other mattress sizes and accessories may be available in the future. For latest info please visit [www.kanmed.se](http://www.kanmed.se)



Kanmed WarmCloud has been used by more than 150 000 patients.

Visit [www.kanmed.se](http://www.kanmed.se) for more information and to find your local supplier for a trial

# KANMED<sup>o</sup>

Kanmed AB  
Gårdsfogdevägen 18B, SE-168 67 Bromma, Sweden  
Telephone: +46 8 564 80 630. Telefax: +46 8 564 80 639  
E-mail: [info@kanmed.se](mailto:info@kanmed.se). [www.kanmed.se](http://www.kanmed.se)